

Promotion!!! 25%OFF!!! Package at \$2,399



German Level Comfort Dental Chair

+



18L Class N Autoclave



Air Compressor



Intraoral Camera



Ultrasonic Scaler



Light Cure



Mini Apex Locator



W&H Handpiece Kit



Cabinet



100pcs gown

QL2028I

Fashionable appearance stable performance

QL2028I latest upgrade model / QL series best-selling model

Fashionable appearance / stable performance

The manufacturing process is further improved / increase the user-friendly functions



Multi-function main control panel

Easy to operate

LED sensor dental lamp

Handy for different
treatment scenarios

Turn the case

Increase the activity space
of assistants

Rotating glass spittoon

The spitting angle is closer
and more convenient

Rotating handpiece hanger

Provide a more convenient operating
angle

External saliva suction device

Convenient for daily filter cleaning

Low-height of chair design

Minimum height **40CM**, which is convenient for kids and the aged, comfortable for dentist to operate.



Flexible rotating case


The rotation angle is 90 degrees, the assistant position has more space for movement, and the four-handed operation is more efficient; It can also facilitate the cleaning and maintenance of the equipment.



Backrest design more suitable for operation

Curved headrest

- The backrest shape design with a wide top and a narrow bottom is more in line with the operating habits of Chinese doctors at 9-11 o'clock. It can approach the patient freely in a more comfortable and healthy sitting posture, and it is not easy to fatigue.
- The curved design provides comfortable support for the neck and prevents the head from shaking left and right.

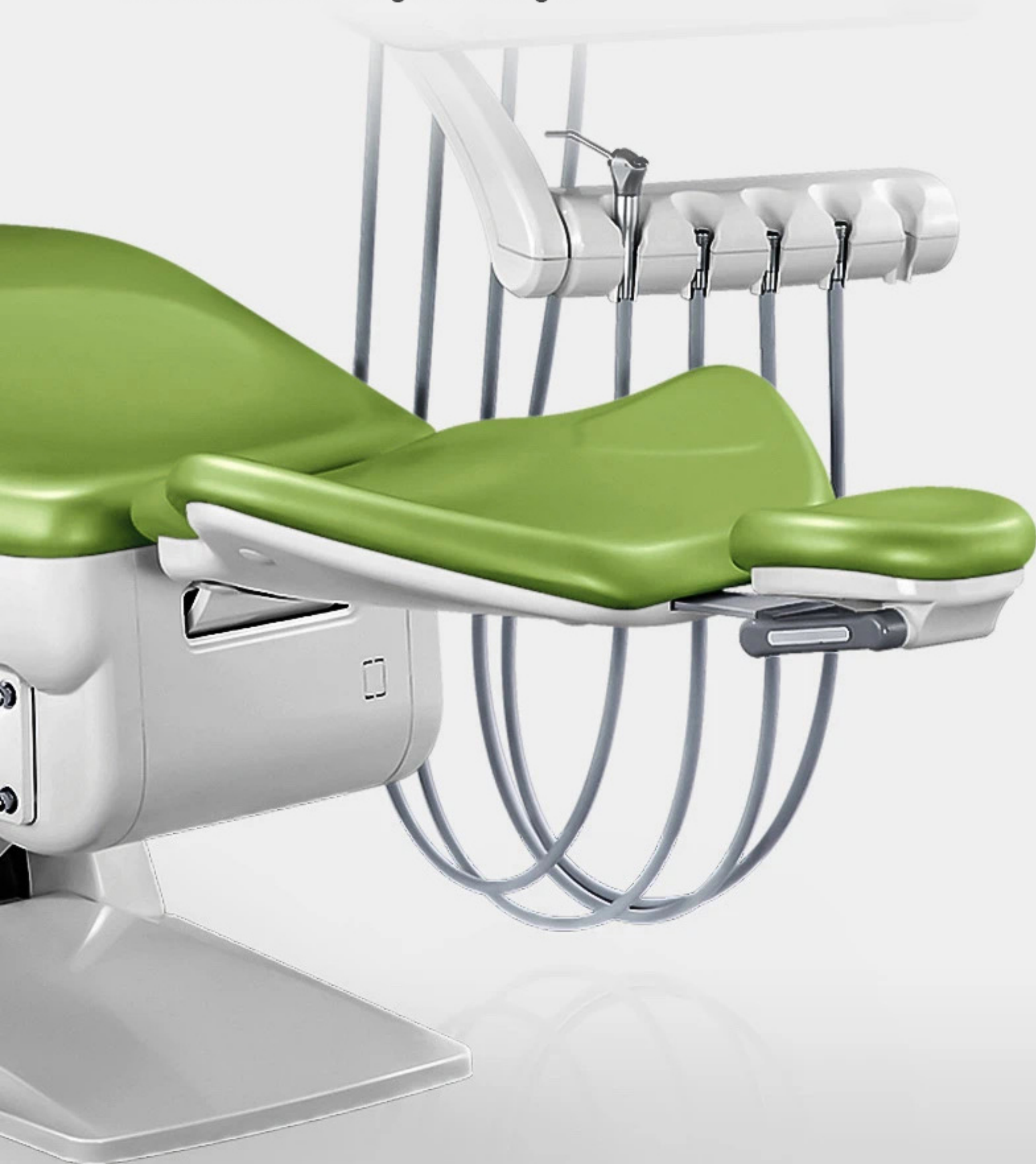

Rotary handpiece holder

- Large size instrument tray, with Silicone pad.
- Rotary handpiece holder can provide convenient angle for operation and avoid direct contact between dentist and handpiece




Adjustable headrest

Headrest can be adjusted ,can quickly adjust position with one hand,
Convenient and efficient



Built out and detachable suction filter

- Screw type tight seal structure, rotating into the lock, medical staff easier to operate
- Suction filter install outside, could see and clean easy.

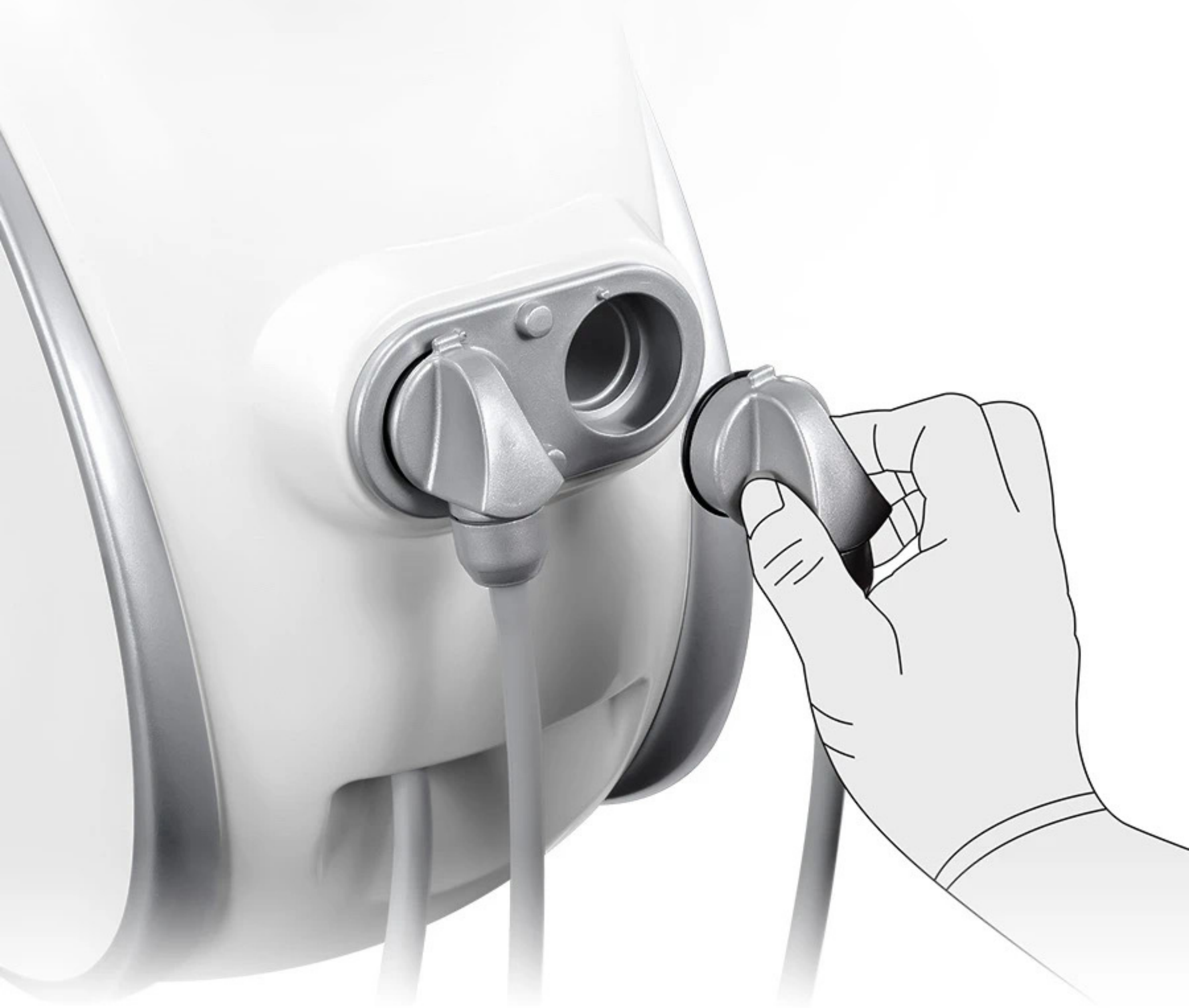

13-key multi-function main control panel

Spitting linkage shortcut key

- Rich functions and simple operation. The postoperative reset function key allows the patient to go up and down more quickly, One key -5 degree emergency position, always guarding safety.
- When the patient needs to spit, the backrest is raised with one button, the spittoon is flushed, and the dental lamp is turned off; After the patient finishes vomiting, the backrest can return to the treatment position with one button, turn off the sputum flushing water, and automatically turn on the dental light; 6 operations are directly simplified into 2 operations, which is more efficient; (the assistant's position is also equipped with a shortcut key for vomiting linkage)



- ①. Setting
- ②. Cuspidor flush
- ③. Cup filling
- ④. Water heating
- ⑤. Zero position
- ⑥. Emergency position (-5°)
- ⑦. Spitting position
- ⑧. Operation light
- ⑨. X-ray film viewer

Rotary glass cuspidor

- Rotary **180°** and detachable spittoon, easy to clean. (Optional spittoon in ceramic)




Combined movable lamp arm

- Aluminum alloy material, light and flexible, stable positioning
- The elastic force could adjust, convenient for upgrade shadowless lamp in the future, Longer service life



Independent pure water supply system

Supply pure water for 3-way spray,high speed handpiece and low speed handpiece of main control system.

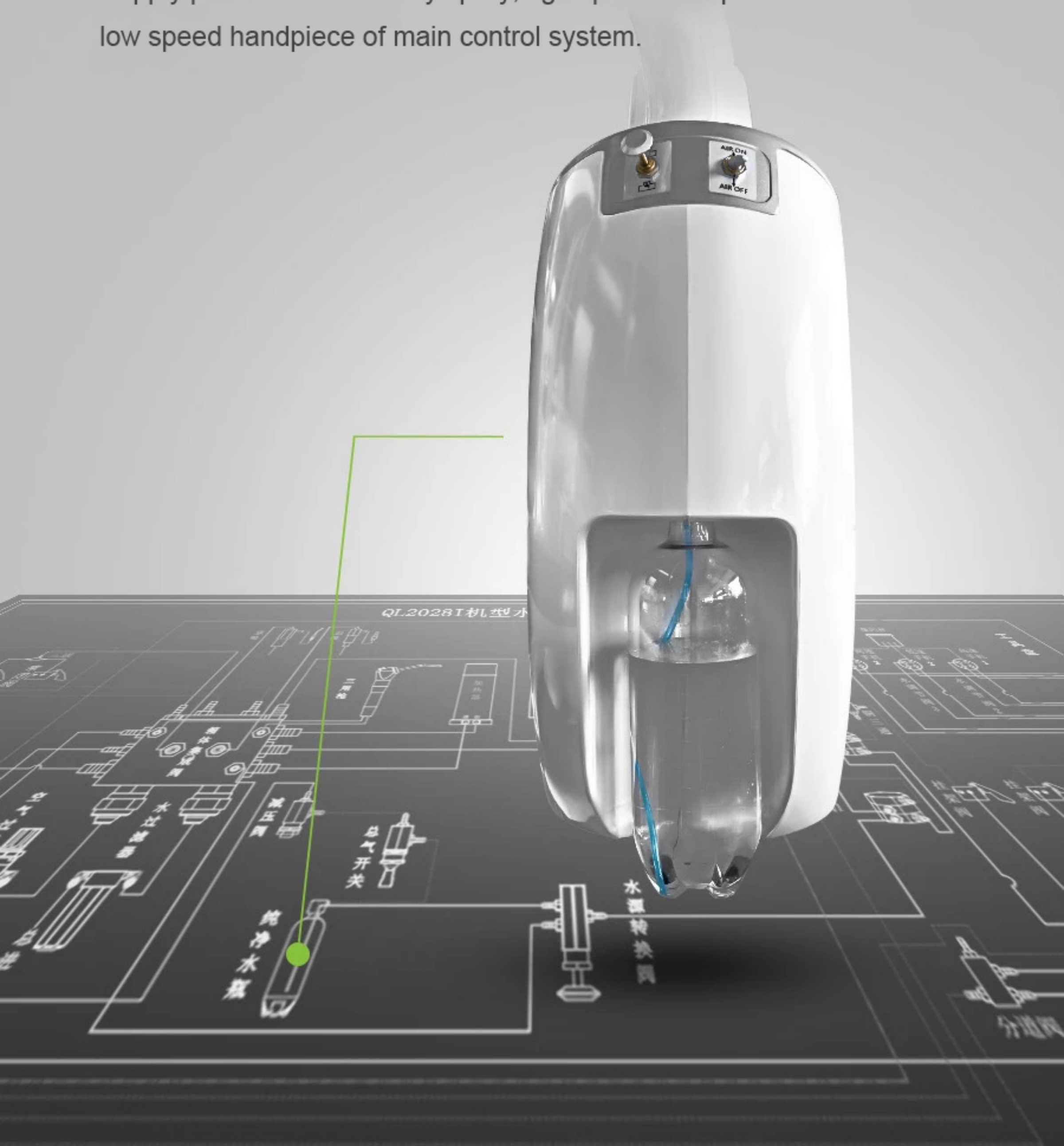

Multi function foot pedal with latest design

Chair position control, with bowl rinsing and cup filling control.



LED sensor light

Two color change,no dazzling

High quality LED illuminating source,
the light spot is soft and uniform,yellow/white/mixed light mode.
Reduce eye irritation to patients, patients are more comfortable
under long-term treatment



Six LED light sources

The light spot is soft, clear and uniform



Sensing and manual control

Take into account different
operating habits



Three modes of luminescence

Yellow/white/mixed light source

