

QL2028 I

牙科综合治疗机 / DENTAL UNIT

轻巧，灵动

Smooth and Streamlined



感应 LED 柔和双色灯

黄光、白光、柔和混合光三种光源模式自由切换，光斑柔和清晰均匀，感应手动两用开关，亮度强至弱可调节，手柄可拆卸，光源分为两组黄光以及四组白光。

高档旋转式玻璃痰盂

痰盂质感强，吐痰舒适，易拆卸清洗，有效防止水流飞溅。

全新下挂器械盘

全新设计，器械盘面加宽，可拆卸便于清洁手机挂架盒，拉手舒适，便于各种操作。

功能性的控制面板

观片灯内置设计，吐痰联动椅位快捷键，口腔灯快捷键，观片灯快捷键，术后复位快捷键，急救 -5°椅位。

吐痰联动椅位

按吐痰椅位快捷键，靠背上升至吐痰位，冲痰功能自动开启，口腔灯自动熄灭，再按此快捷键，则自动恢复刚才工作状态的椅位 冲痰关闭，口腔灯开启。

全新设计多功能脚踏

独特的金属支架，重新定义万向开关，根据人体工程学设计，使用更舒适，具有椅位冲痰供水控制功能。

LED sensor light with yellow and white light mix

More tender light: 6 light source with 4 white and 2 yellow; Easy switch of yellow light mode, white light mode and mix light mode; Super clear light spot, brightness adjustable from strong to weak; Detachable and sterilizable handle.

Rotary glass cuspidor

Rotary and detachable cuspidor, easy to clean.

Hanging delivery tray

New design, detachable & easy clean handpiece holder block, and comfortable tray handle.

Full function control panel

X-ray film viewer built in design, with spitting position setting, operation light control, x-ray film viewer control, zero position control, Emergency position(-5°) control.

Spitting position setting

Press the spitting position button, automatically backrest upward, bowl rinsing start, light turns off; press again, automatically return to working status, backrest backward, bowl rinsing closes, light turns on.

Multi function foot pedal with latest design

Unique metal handle, with chair position control, with bowl rinsing and cup filling control.



种植手术口腔灯

LED light

- 全封闭设计，配轻松拆卸的高温消毒的把手套，可清洁和消毒；
- 该灯具有可调色温（4000–5500K，7档），高亮度（可达 60000LUX），无影效果好；光斑方正（85*160）
- 是一款既可以用于种植手术照明又可以兼顾其他口腔手术照明多用途口腔灯。
- Totally enclosed design and hygienic autoclavable handle covers which is easily taken out or put in
- 7color temperatures to be chosen(4000K to 5500K) and high illumination up to 60000Lux.
- The light is also an excellent choice for the ENT and cosmetic surgery.
- (rectangle pattern,high illumination, superb shadow less effect)



医生椅 / DOCTOR CHAIR



A6



A5



A4



A3

脚踏

MULTIFUNCTIONAL FOOT PEDAL

全新设计多功能脚踏
Multi function foot pedal with latest design

独特的金属支架，重新定义万向开关，根据人体工程学设计，使用更舒适，具有椅位冲痰供水控制功能。
Unique metal handle, with chair position control, with bowl rinsing and cup filling control.



可选配置

OPTIONS



工具盘 / Operation tray



脚踏 / Multifunctional foot pedal

选配部件 / OPTIONS



台湾堤摩讯电机 / Taiwan Timotion motor



光固化 / Curing light



洁牙机 / Scaler



内窥镜系统 / Endoscope system



陶瓷痰盂 / Ceramic cuspidor



内窥镜支架 / Screen holder


无缝座垫 / PU LEATHER



超细纤维皮座垫 / MICROFIBER LEATHER



安装示意图 / INSTALLATION DIAGRAM

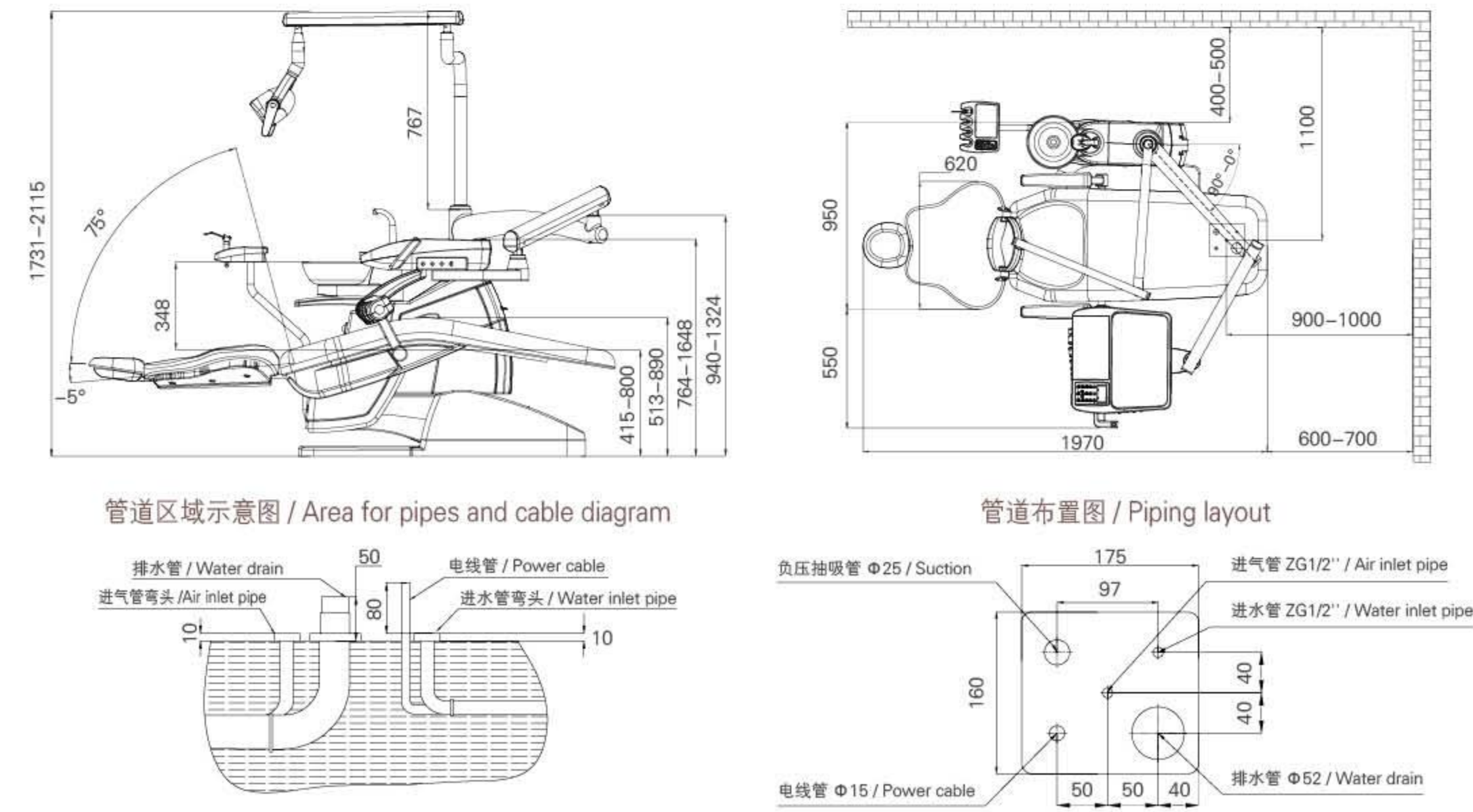

牙椅配置表 / CONFIGURATION OF DENTAL UNIT

类 目 Classification	配 置 Configuration	型号 Model				
		QL2028 IV	QL2028 III	QL2028 II	QL2028 I	QL2028
治疗椅 Patient Chair	豪华的超纤皮座垫	★	☆	☆	☆	☆
	无缝 PU 皮座垫	☆	★	★	★	★
	联动椅位补偿功能	★	★	★	☆	☆
	机椅互锁系统	★	★	★	☆	☆
	水气电一键开关装置	★	★	☆	☆	☆
	椅位记忆功能（9 组）	★	★	☆	☆	☆
	吐痰复位	★	★	★	★	★
	急救椅位	☆	☆	★	★	★
	豪华型多功能脚踏	★	★	★	★	★
口腔灯 Operation Light	豪华型 LED 冷光灯	★	★	★	★	★
器械盘 Operation tray	上挂式	★	☆	☆	☆	☆
	下挂式	☆	★	★	★	★
	十六键多功能控制面板	★	★	☆	☆	☆
	十三键多功能控制面板	☆	☆	★	★	★
	外置观片灯	★	★	☆	☆	☆
	内置观片灯	☆	☆	★	★	★
	三用枪（冷）	★	★	★	★	★
主机箱 Main side box	可旋转拆卸玻璃痰盂	★	★	★	★	☆
	强吸和弱吸	★	★	★	★	★
	自动冲痰供水恒温控制系统	★	★	★	★	★
	纯净水供应系统	★	★	★	★	★
副控台 Assistant control	八键多功能操作面板	★	★	★	★	★
	三用枪（热）	★	★	★	★	★
	可调节强吸接头	★	★	★	☆	☆
地箱 Floor box	内置型	★	★	★	☆	☆
	外置型	☆	☆	☆	★	★
医生座椅 Doctor's stool	豪华型	A5	A5	☆	☆	☆
	标准型	☆	☆	A4	A3	A3
选配部件 Options	高速手机	▲	▲	▲	▲	▲
	低速手机	▲	▲	▲	▲	▲
	内窥镜系统	▲	▲	▲	▲	▲
	超声波洁牙机	▲	▲	▲	▲	▲
	LED 光固化机	▲	▲	▲	▲	▲
	医用空气压缩机	▲	▲	▲	▲	▲
	抽吸机	▲	▲	▲	▲	▲
	上挂式工具盘	▲	▲	▲	▲	▲
	托盘水杯筒纸巾盒	▲	▲	▲	▲	▲
	豪华超纤皮座垫	▲	▲	▲	▲	▲
	进口电磁阀	▲	▲	▲	▲	▲
	台湾堤摩讯电机	▲	▲	▲	▲	▲

备注 NOTE: ★ 标准配置 Equipped / ☆ 无此配置 None / ▲ 选配部件 Options