

Configuration	★ Standard / △ Optional	Luxury	Standard				Luxury	Standard
Patient Chair	Luxury microfiber leather cushion	★	★	Main Side Box	90 degrees rotatable side box	★	★	
	Soft-start and soft-stop function & backrest anti-collision	★	★		Rotatable Ceramic Cupidor	★	★	
	Synchronized Chair & Chair Interlock System	★	★		Imported water & air tube, two water bottles	★	★	
	Spitting / Treatment / Emergency Chair Position	★	★		Independent water filter, dual filtration function	★	★	
	Luxury Multi-functional Foot Pedal	★	★	Control System	HD LCD screen, air pressure and rotating speed display	★	★	
Instrument Tray	Instrument Tray heating Water System	★	★		Built-in LCD X-ray film viewer	★	★	
	6-position Adjustable Instrument Holder	★	★		9 sets of chair position memory	★	★	
	One button Water-air-electricity on/off switch	★	★		Intelligent self-test, fault code display	★	★	
	Handpiece independent water adjustment	★	★	Pipeline Disinfection	Built-in disinfection module	★	★	
	Turnable stainless steel tray	★	★		One-click automatic disinfection system / 9 minutes	△	★	
Assustant Control	3-way Syringe (with warm water)	★	★		One-click automatic disinfection system / 25 minutes	★	△	
	9 keys control panel	★	★	Water volume test, one-button water source switching	★	△		
	Adjustable Strong/Weak Suction Head	★	★	Options	Selected imported leather cushion	△	△	
Floor Box	Built-in Floor Box	★	★		Built-out floor box	△	△	
	Voltage stabilizing filter and Voltage display	★	★		High and low speed handpiece	△	△	
Dentist Stool	Comfortable Assistant stool	★	★		Built-in scaler and Built-in curing light	△	△	
	Professional Dentist stool	★	★		Intra-oral camera and its holder	△	△	
Dental Light	8 bulbs LED sensor light (3 kinds of light mode)	★	★		Air Compressor	△	△	

M200



HIGH-END
DENTAL CHAIR

INNOVATIVE IDEA, EXQUISITE CRAFTSMANSHIP

MIPONT Medical is located in Nansha District, a key development area of Guangzhou City, with vast development potential and a solid industrial foundation. The company adheres to strict product inspection and process standards, with advanced production and quality control systems. As an emerging enterprise in dental chair manufacturing, MIPONT Medical is committed to innovation and development, providing high-quality and refined solutions for dental clinics.



DESIGNED FOR MODERN HIGH-END DENTAL CLINICS




TO CREATE A MODERN HIGH-END CLINIC

Has successfully develop dental unit tailored to the market, with certification of innovative enterprise and intellectual property rights, relying on technological innovation, obtaining competitive advantages in the market for sustainable development.

Designed for modern high-end dental clinics

Adhering to the design philosophy of "Less is More", MIPONT emphasizes the combination of elegant design and medical needs, infusing products with innovative designs tailored for modern high-end dental clinics, with strong insight of user needs and market needs.

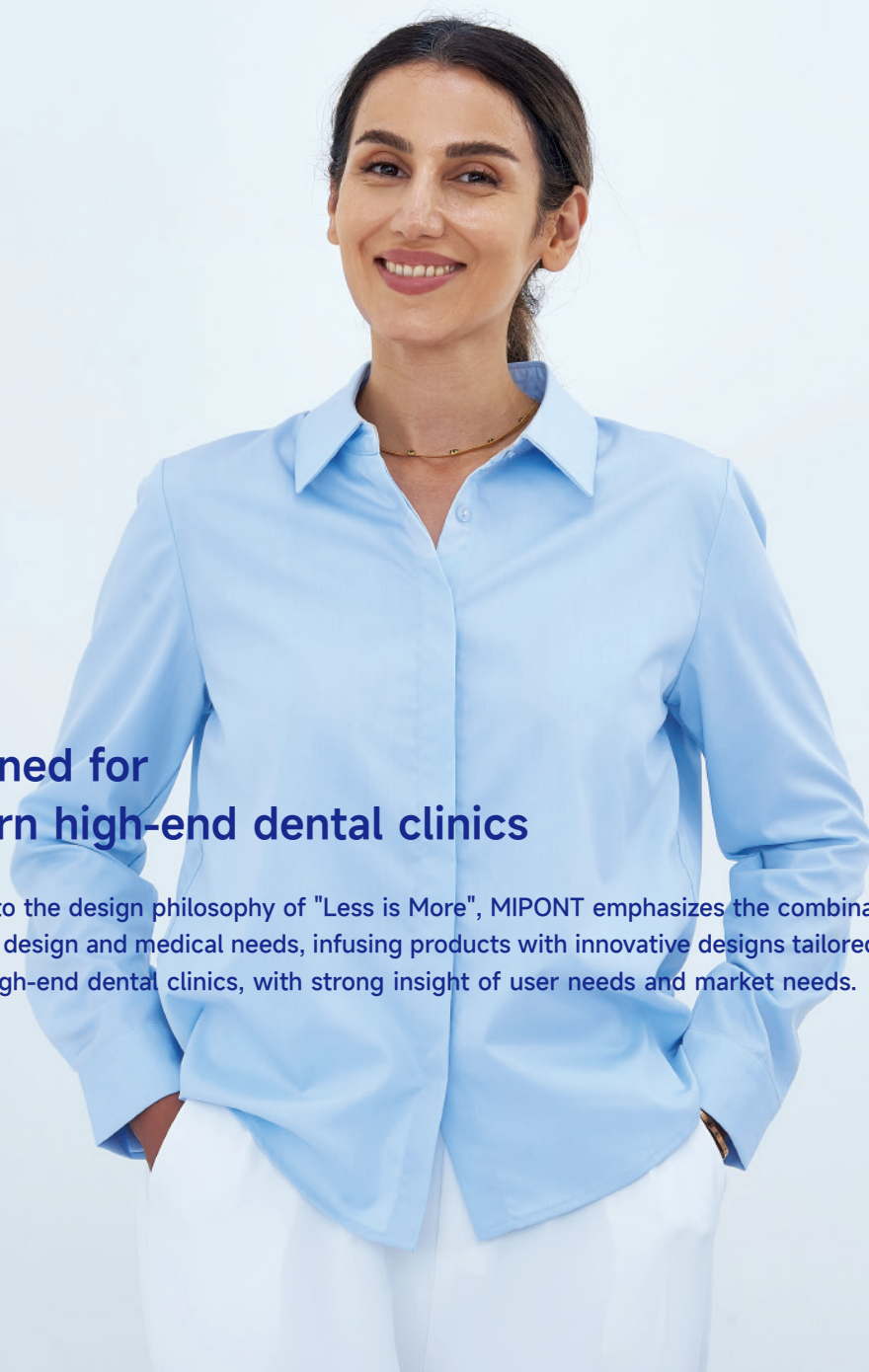

Has established professional lab with a variety of lab equipment and over 30 sets of lab equipment, running over 150 tests covering more than 280 kinds of materials, including environmental testing and reliability testing and etc.

- OVER 30 Professional lab equipment
- OVER 50 Professional Quality Control team
- OVER 150 Testing standards

Production and Management

Has a strict and lean production management system and a refined 5S management workshop to ensure the efficient operation of production, providing high quality product and dedicated service to our partners.

Strict Quality Control

Conducts steps of quality control from the raw material with product reliability test, simulation test of the use of dental chairs, extensive research and testing of the existing products, and product upgrades, to ensure well-designed and well-produced high-end product.





M200 Standard

Smart and high-end, empowering precise treatment.

M200 Standard is designed with smart control system and leading configurations, which is perfect solution tailored to meet the real needs of doctors, patients, and assistants

Experience the unparalleled comfort.



PROVIDES SUPERIOR COMFORT EXPERIENCE

High-elastic microfiber leather cushion, ensuring long-lasting comfort.

380mm chair position, ideal for elderly and pediatric patients.
Soft-start and soft-stop function for a smooth experience.

SOFT AND FINE

Microfiber leather cushion, increase patient comfort.
Front-flip armrests with leather-covered padding, delicate touching.

▼ Standard 2 stools configuration

MIPONT's in-depth study of ergonomic design, comfortable lumbar support.
The stool's height and backrest recline are fully adjustable for optimal comfort.

Comfortable assistant's stool
Luxury dentist's stool



Luxury dentist's stool with microfiber leather and lumbar support.
Comfortable assistant's stool with quiet wheels for smooth mobility.

Superior configuration.



MULTIDUNCTIONAL FOOT PEDAL

7 integrated functions, effortless and efficient operation.
Delivers a new level of precision.

Independent water and air control to prevent cross-infection.
Spittoon linkage simplifies operation with a single touch.

OUTSTANDING LIGHTING



LED HIGH-END SENSOR DENTAL LIGHT

Soft lighting with 8 imported PHILIPS LEDs, reducing direct eye irritation for patients.
Two operating modes to meet the needs of different treatment scenarios for dentists.
100,000 high-frequency life cycle tests ensure quality, safety, and durability.

Manual / sensor dual-mode precision control, with adjustable brightness ranging from 8 to 40k Lux.
Three light source modes for easy switching: yellow light, white light, and yellow-white mixed light.

Pro smart disinfection system

STRICT CONTROL FOR DAILY DISINFECTION

Simplified process with one-click disinfection operation.

visualized precise control of the disinfection progress.
Full pipeline coverage for the water system.



Tap water rinse



Disinfectant rinse



Disinfectant retention



Tap water rinse



Purified water rinse



MOS Operation System

KNOW YOU BETTER

Accurate and easily control of various functions

10 self-check function at a glance, display the fault code if any
High-definition LCD screen, digital air pressure and speed display
9 groups of chair positions, providing custom chair position settings

System self-check

Assistant control board communication
Main control board communication
Foot control board communication

Safety lock function

Dental light function

X-ray viewer function

Water heater function

Water supply function

Spittoon flushing function

Obstruction detection function





M200 Luxury

Smart and intelligent disinfection.

MIPONT M200 Luxury is a dental unit with one-key intelligent disinfection function, running automatically without more operation and covering all pipes disinfection.

“ Simply the best of
intelligent disinfection dental unit ”



MAX

ONE-KEY INTELLIGENT DISINFECTION SYSTEM

Double-bottle design

Water can be easily filled from the top of the lid without opening the bottle.

Built-in disinfection module

Built-in module in the sidebox, easy to use. Main pipeline covered, fully disinfect the entire pipeline and the strong & week suction.

Intelligent emergency response function

After the automatic disinfection is completed, the unit will automatically turn off the water, electricity and gas without any operation.

One-click intelligent disinfection

Automatically disinfection system without supervision. Disinfection interruption recovery, power failure recovery and water shortage reminder.


2 stools configuration

Comfortable nurse stool

Silent pulley, smooth movement.

Luxury dentist stool


Microfiber leather provides comfortable lumbar support.

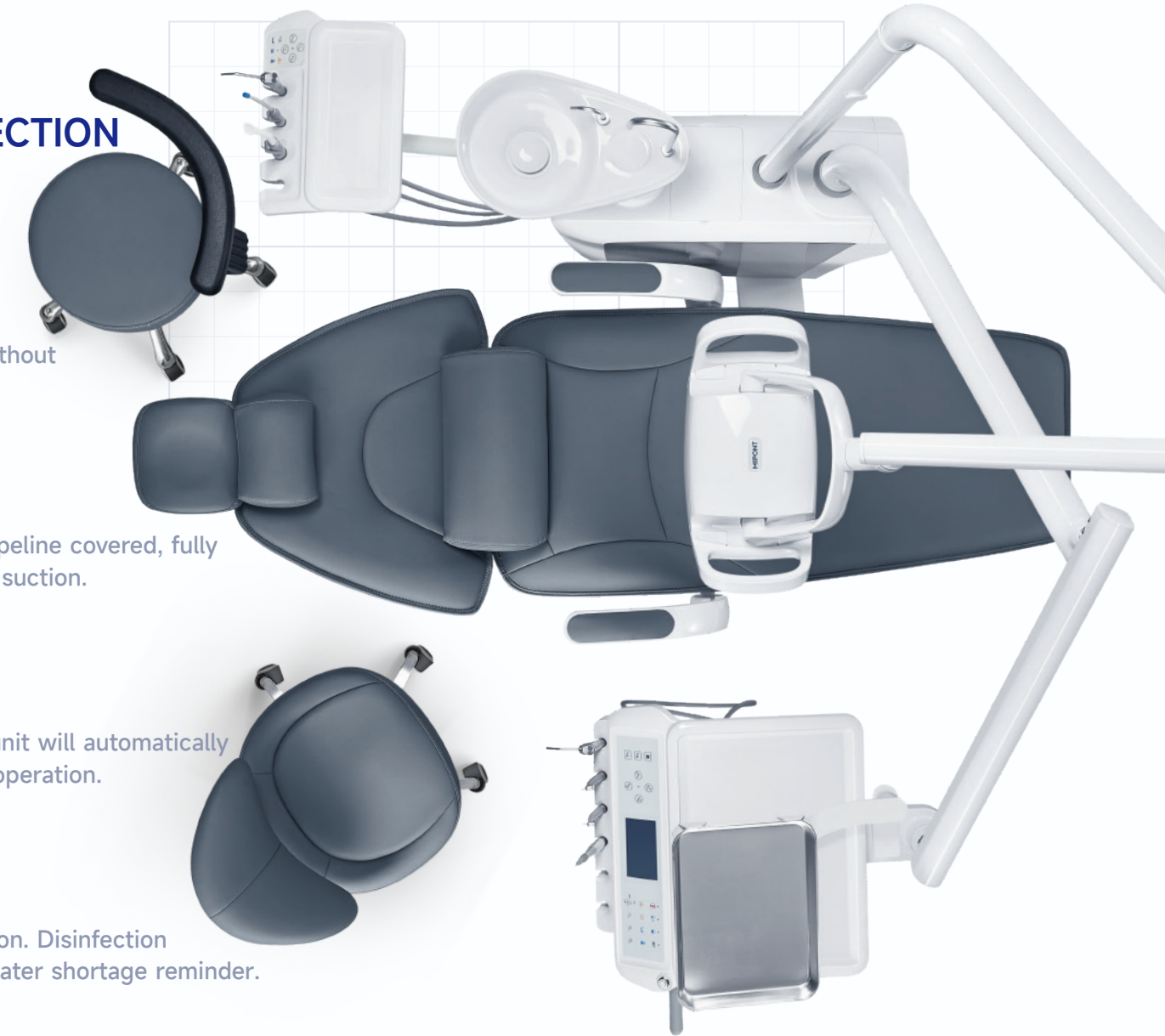

MOS

Operation System

M.

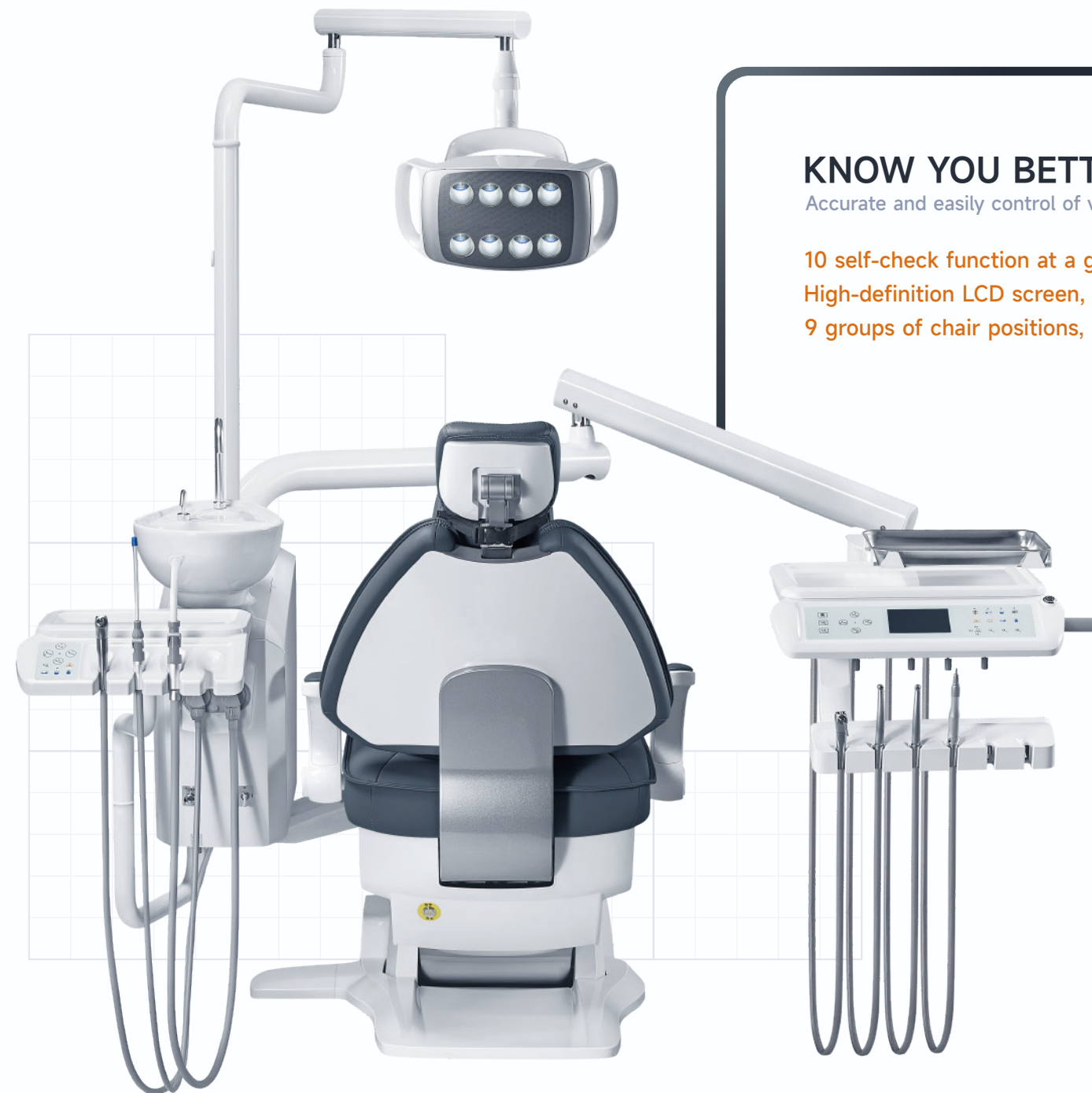

KNOW YOU BETTER

Accurate and easily control of various functions

10 self-check function at a glance, display the fault code if any
High-definition LCD screen, digital air pressure and speed display
9 groups of chair positions, providing custom chair position settings

System self-check

Assistant control board communication

Main control board communication

Foot control board communication

Safety lock function

Dental light function

X-ray viewer function

Water heater function

Water supply function

Spittoon flushing function

Obstruction detection function

Long-term Solution For Clinic

Additional hot water function in instrument tray

Provide warm water for handpiece and syringe, effectively reducing pain.

Safety with anti-collision system

Emergency stop makes the operation safer when encountering obstacles.

Imported motor

Leading quality, quiet and stable.

Imported water and air tube

Good flexibility, corrosion resisting, durable and high quality.

Additional tray

Stainless steel is more durable, corrosion-resistant and sterilizable.

Digital voltmeter

Display of input voltage number, additional filter, anti-voltage fluctuation.



EXCELLENT DENTAL LIGHT, BRAND NEW EXPERIENCE

HIGH-END LED DENTAL LIGHT

Gentle light source reduces direct stimulation to patients' eyes.
100,000 times high-frequency life test.
8 imported lamp beads from PHILIPS.

Manual/sensing dual-mode for different treatment

3 modes: yellow, white, and mixed mode

8-40k Lux adjustable brightness

Strong metal backrest

3MM thickened full metal base for stronger support.

Soft Start and Soft Stop

Chair movement without sudden quiver.
More comfortable for patient.

Thicker curved plate, firmer connection

Extra large die-cast 30mm aluminum plate.

Solid bottom plate

Thickness 12mm, stable as a rock.
The unit weighs 250KG, making it more stable.



Center connection to the sidebox

Easier to install the equipment outward rotatable sidebox
Provide flexibility for narrow space installation without splitting the sidebox.



Different Color
Different Style

Microfiber (Stand for M200)

					
M06 Dark Blue	M03 Blue	M07 Grass Green	M05 Apple green	M01 Beige	M08 Yellow
					
M02 Red	M10 Dark Red	M09 Purple	MC11 Dark Grey	MC12 Light Grey	M04 Brown

Specially Selected Microfiber Leather (Optional)

					
MF01 Bluish Green	MF02 Apple Green	MF03 Blue	MF04 Dark Blue	MF05 Rose Red	MF06 Dark Red