

Configuration

V3 Luxury dental unit
Reamrks: ★standard / ☆optional parts

Main side box

- ★ One key automatic disinfection system
- ★ Rotary side box
- ★ Rotary and removable glass cuspider
- ★ water cup induction automatic water supply
- ★ Imported water & air tube
- ★ Two water storing bottles
- ★ Independent water filter, dual filtration

Instrument tray

- ★ 3 doctors with 9 sets of chair position memory
- ★ 6-position Adjustable instrument holder
- ★ Handpiece Anti-retraction function
- ★ Water-air-electricity with one on/off switch control
- ★ Handpiece independent water adjustment function
- ★ Built-in X-ray film viewer
- ★ Small tray
- ★ Original anti-fouling silicone pad

Patient chair

- ★ Skin-friendly leather cushion (high-quality silicone)
- ☆ Premium PU leather cushion
- ★ Soft Start and Soft Stop function
- ★ Synchronized chair design & Chair interlock system
- ★ Spitting chair position
- ★ treatment chair position
- ★ emergency chair position
- ★ Backrest anti-collision
- ★ Taiwan Timotion motor
- ★ Luxury multi-functional foot pedal

Operation light

- ★ Direct operation light (Imported 6 lamp beads, with induction)
- ☆ Reflected operation light (yellow and white light, with induction)

Assistant tray

- ★ 3-way syringe (hot)
- ★ 9 keys control panel
- ★ Adjustable strong/weak suction
- ★ Original anti-fouling silicone pad

Dentist stool

- ★ Luxury dentist stool (skin-friendly leather)
- ☆ Professional dentist stool (High class PU leather)
- ★ Luxury nurse stool (skin-friendly leather)
- ☆ Professional nurse stool (High class PU leather)

Floor box

- ★ Built-in floor box
- ★ Digital voltage display
- ★ Equipped with wave filter

Options

- ☆ Upgrade to mobile cart type
- ☆ Built-out floor box
- ☆ High and low speed handpiece
- ☆ Built-in scaler
- ☆ Built-in curing light
- ☆ Endoscope system, Endoscope holder
- ☆ Air Compressor

Color Chart

V3 Luxury dental unit
Skin-friendly Cushion

V01 sky blue

V02 light blue

V03 medium blue

V04 light gray

V05 greenish brown

V06 light red

V07 light orange

V08 apple green

TO BUILD A CLINIC WITH HIGH LEVEL

V3
dental unit
V3 Luxury



WEILAN

V3

TO BUILD A CLINIC WITH HIGH LEVEL

Professional dental equipment manufacturing
enterprises to build high-end dental chairs



V3

V3 LUXURY




EASY TO DISINFECT

Makes the disinfection

EASIER AND MORE EFFICIENT

Cover all tubes

Disinfect to all water tubes, worry-free sensory control.
It can realize the disinfection of handpiece, 3 way syringe, cup filling tube, and strong & weak suction tubes.

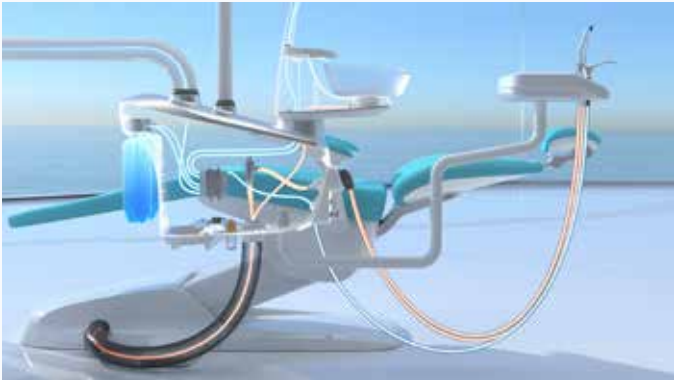

One key automatic disinfection

Completing the disinfection operation with one button,getting rid of the traditional complicated type.

Automatic operation is faster

Modular system, disinfection procedures intelligent operation.
The disinfection time is accurately controlled to ensure that the disinfection effect reaches to standard each time.

LEADING TECHNOLOGY

Induction water supply

Automatic sensing controls the water supply.
Gravity sensing precisely controls the water volume.



Exclusive Design Dental Light

- 6 imported lamp bulbs
- 3 light source modes

Simple technology, sensor and button dual control mode.
White, yellow, and mixed light source to meet the needs of different treatment operations.



Self-developed multifunctional pedal

Comfort control

Optimize the control angle and feedback strength, comfortable control, relaxing work.

Fully upgraded functions

Easy to control the dental light, the reset of treatment position, spitting position, and emergency chair position, improving treatment efficiency.

More efficient one-key dual function

- Combine the spitting and reset position control as one
- Combine the dental light and chair position control as one

Improve treatment efficiency

HIGHER CONFIGURATION

Higher class-leading configuration


High quality materials to enhance quality

The whole dental chair adopts SMC medical-grade imported water-air pipe, makes it more durable.
Double water filtration, improve water quality, protect the valve; double pure water bottles, water source free switch, larger capacity.

More convenient four-hand operation

Equipped with professional doctor's chair and nurse's chair.
Rotatable main box, providing more operating space for nurses.

- Nurse's chair
Silent pulley, high-density shaping cotton, comfortable sitting.
- Doctor's chair
Backrest waist protection design, fit to support the waist, optimize the curve of the cushion, relieve the pressure of the legs.

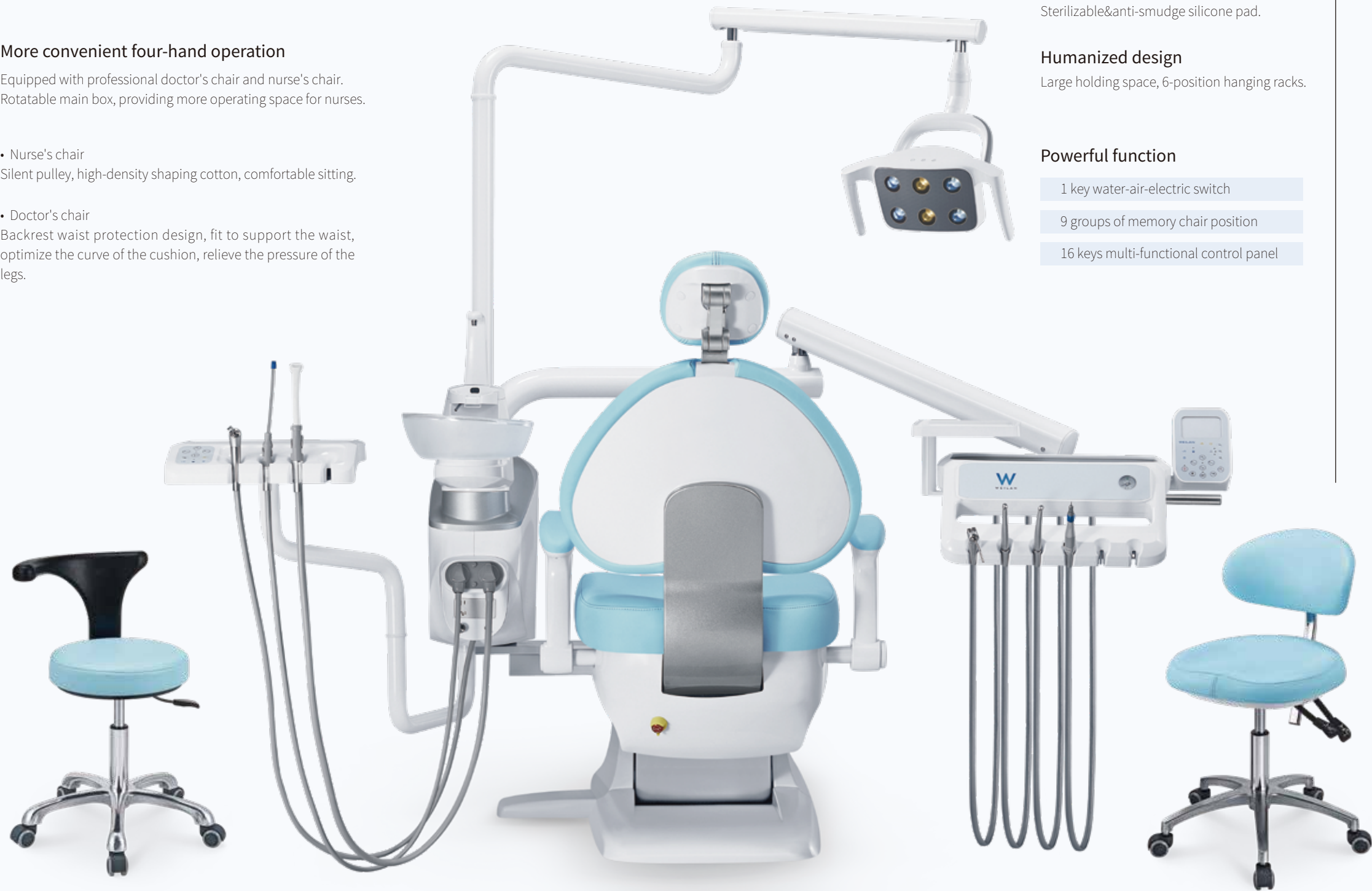

Instrument tray

Focus on infection control

Handpiece anti-suckback design.
Sterilizable&anti-smudge silicone pad.

Humanized design

Large holding space, 6-position hanging racks.

Powerful function

- 1 key water-air-electric switch
- 9 groups of memory chair position
- 16 keys multi-functional control panel

COMFORT CREATION



Ultimate skin-friendly experience

Skin-friendly Cushion

Skin-friendly touch, ultimate skin-friendly experience.
High-quality memory foam with automatic shaping, comfortable lumbar support neck.

Comfortable backrest

One-piece metal backrest brings calm and support.
No obstruction to the doctor's legs, free access to the patient.

Curved headrest

Curvilinear curvature to fit the head.
Adjustable multi-joint angle for one-handed operation and easy positioning.

Soft start / stop

Lightweight soft start-stop, significantly reducing the stuttering and back-rubbing sensation during the start-stop process, enhancing the patient's comfortable experience during the operation.

V3 IMPLANT



Bring the implant operation with more concentration and ease



Exclusive Design Implant Light

- 21 high CRI imported light sources bulbs
- 3 gears general color temperature/light color (K): 4200K / 4700K / 5200K
- One key anti-curing mode
- CRI>90Ra
- Light life >40000h

- Implant surgery mode

Target distance 0.7 meters Illuminance value (Lux)
10 levels of illumination adjustable: 8000-80000lx

- Dental diagnosis mode

Target distance 0.7 meters Illuminance value (Lux)
10 levels of illumination adjustable: 6000-25000lx

Aluminum alloy detachable handle, easy to disinfect.
Separate control panel, induction control switch, no contact sensory control.
Excellent lighting brightness, adjustable in ten levels, any surgical area in the oral cavity can be seen clearly.
21 imported light source, double light spot mode, providing excellent shadowless lighting effect at any angle.



Mobile trolley instrument tray

Electronically controlled lift, adjustable and suitable for any treatment height. 220v independent power supply, meeting the power requirements of various treatment equipment.